

NOTES

Rapid Detection of Vaginal *Candida* Species by Newly Developed Immunochromatography^{∇†}

Hidehito Matsui,¹ Hideaki Hanaki,^{1*} Kengo Takahashi,² Akihiko Yokoyama,² Taiji Nakae,¹ Keisuke Sunakawa,³ and Satoshi Omura^{1,3}

Laboratory for Antimicrobial Agents, Kitasato University, 5-9-1 Shirokane, Minato-Ku, Tokyo 108-8641, Japan¹;
 Isehara Research Laboratory, Kanto Chemical Co., Inc., 21 Suzukawa, Isehara, Kanagawa 259-1146, Japan²;
 and Graduate School for Infection Control Sciences, Kitasato University, 5-9-1 Shirokane, Minato-Ku, Tokyo 108-8641, Japan³

Received 13 May 2009/Returned for modification 22 June 2009/Accepted 23 July 2009

For the diagnosis of vulvovaginal candidiasis, we developed a simple immunochromatographic method that enables the detection of vaginal *Candida* spp. within about 30 min. Overall, the sensitivity, specificity, positive predictive value, and negative predictive value of this method appeared to be 80.3, 99.3, 98.0, and 92.0%, respectively.

Vulvovaginitis is one of the most common infectious diseases in women's clinics. In the United States, 6 to 10 million cases of gynecology clinic visits per year were estimated to be due to vaginitis (4, 7). The major etiologies of vaginitis are infection by bacteria, fungi, and trichomonads (10). Among vulvovaginitis cases caused by fungi, about 80 to 90% are due to infection with *Candida albicans* (8, 11).

In routine clinical practice, vulvovaginal infection is primarily diagnosed from the patient's complaints and clinical symptoms such as pruritus and increased vaginal discharge. Symptomatic vulvovaginitis may be readily diagnosed from microscopic examination of the vaginal discharge. However, this method can detect only 40 to 70% of *Candida* species compared with the culture method (5, 10) and requires experience and the use of expensive

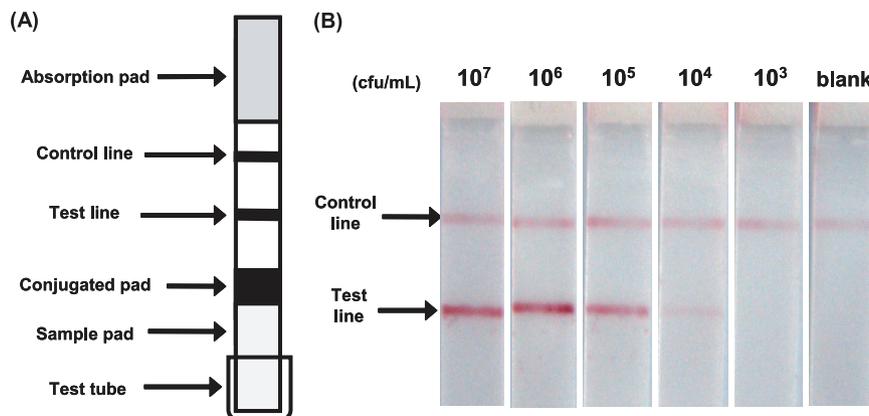


FIG. 1. Immunochromatographic test strip. (A) Schematic representation of the strip. The nitrocellulose membrane strip consists of anti-mannan IgG and anti-rabbit IgG antibodies that are immobilized as a test line and a control line, respectively. The conjugated pad contains the colloidal gold–anti-mannan IgG conjugate. The proximal and distal ends of the nitrocellulose membrane are covered with the sample pad and the absorption pad, respectively. (B) Representative results. Suspensions of *C. albicans* ATCC 10231 at cell densities of 10^3 to 10^7 CFU/ml were heated in antigen extraction buffer and subjected to immunochromatography for 20 min. The sample with double red lines is *Candida* positive, and that with a single red line only at the control line is negative. Data for other *Candida* species are not shown. Only a portion of the strip is shown.

* Corresponding author. Mailing address: Laboratory for Antimicrobial Agents, Kitasato University, 5-9-1 Shirokane, Minato-Ku, Tokyo 108-8641, Japan. Phone: 81-3-5791-6413. Fax: 81-3-5791-6415. E-mail: hanaki@insti.kitasato-u.ac.jp.

† Supplemental material for this article may be found at <http://cvi.asm.org/>.

∇ Published ahead of print on 5 August 2009.

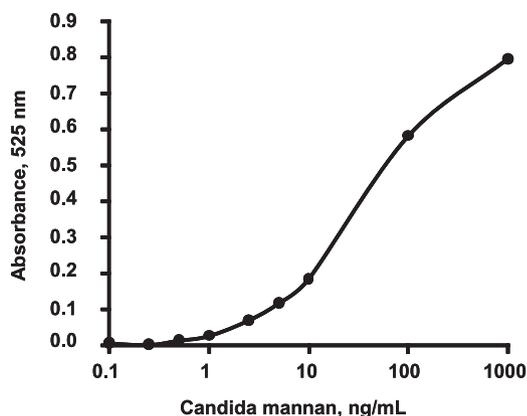


FIG. 2. Correlation between the amount of *Candida* mannan and the intensity of the test line. Purified *C. albicans* mannan was serially diluted and then subjected to immunochromatography. The test line was traced by a densitometer at 525 nm. The absorbance was plotted against the amount of mannan quantified by the phenol sulfuric acid method with D(+)-mannose as the reference (3).

equipment. The culture method may be regarded as the standard; it requires skill, costly equipment, and several days to obtain a reliable result. Accordingly, a rapid and simple method to detect vulvovaginal *Candida* infection has been long awaited. We developed an immunochromatographic method that can detect vulvovaginal *Candida* infection within about 30 min.

The immunochromatographic strip was prepared with anti-*Candida* mannan polyclonal antibody (anti-mannan immunoglobulin G [IgG]) and anti-rabbit IgG antibody serving as the test line and the control line, respectively (Fig. 1A; see the

supplemental material). The colloidal gold-anti-mannan IgG conjugate was placed on a conjugated pad. The immunochromatography works as follows. A test strip (Fig. 1A) is kept standing in a small test tube containing 0.5 ml of the sample to be tested, and the chromatography is developed at ~24°C for 20 min. The sample developed in the sample pad reaches the conjugated pad, and the *Candida* mannan in the solution reacts with the colloidal gold-anti-mannan IgG conjugate. This complex further migrates toward the distal end of the test strip (Fig. 1A). When this complex reaches the test line, immobilized anti-mannan IgG traps the complex and forms a sandwich-type immunocomplex. The product appears as a red line (Fig. 1A, test line). The colloidal gold-anti-mannan IgG conjugate that has not reacted with *Candida* mannan passes the test line and reaches the control line. The anti-rabbit IgG antibody immobilized at the control line traps the conjugate, forming a second red line (Fig. 1A, control line).

First we tested the sensitivity of immunochromatography with purified *Candida* mannan that was extracted from *C. albicans* ATCC 10231 and quantified as reported previously (3, 9). Extracts containing a *Candida* mannan level over the threshold of 0.5 ng/ml showed a positive reaction at the test line. The test line was scanned with a densitometer, and the absorbance was plotted against the amount of *Candida* mannan. The absorbance curve appeared as a pseudolinear line in a concentration-dependent manner (Fig. 2). *Candida* mannan at concentrations below about 0.5 ng/ml was undetectable macroscopically.

In the next experiment, the *Candida* suspension was adjusted to a cell density of 10³ to 10⁷ CFU/ml and *Candida* mannan was extracted by heating the cells in antigen extraction

TABLE 1. Strains used in this study and their culture conditions

Organism(s)	Culture conditions
<i>Candida</i> species: <i>Candida albicans</i> ATCC 10231, <i>Candida tropicalis</i> ATCC 750, <i>Candida guilliermondii</i> ATCC 6260, <i>Candida parapsilosis</i> ATCC 22019, <i>Candida glabrata</i> ATCC 2001, <i>Candida krusei</i> ATCC 6258.....	Sabouraud dextrose agar, 35°C, 48 h
Other fungal species: <i>Trichosporon cutaneum</i> IFM40140, <i>Trichosporon asahii</i> IFM48575, <i>Trichosporon montevidense</i> IFM51966, <i>Trichosporon mucoides</i> IFM48611, <i>Cryptococcus neoformans</i> clinical isolate no. 365, <i>Cryptococcus curvatus</i> clinical isolate no. 334, <i>Rhodotorula mucilaginosa</i> IFM48529, <i>Saccharomyces cerevisiae</i> IFM40022.....	Sabouraud dextrose agar, 35°C, 48 h
Bacterial species	
<i>Staphylococcus aureus</i> FDA209P, <i>Staphylococcus epidermidis</i> clinical isolate no. 1, <i>Micrococcus luteus</i> ATCC 9341, <i>Enterococcus faecalis</i> ATCC 29212, <i>Enterococcus faecium</i> NCTC12204, <i>Escherichia coli</i> NIHJ JC-2, <i>Citrobacter freundii</i> ATCC 8090, <i>Klebsiella pneumoniae</i> NCTN9632, <i>Klebsiella oxytoca</i> clinical isolate no. 1, <i>Proteus mirabilis</i> IFO3849, <i>Proteus vulgaris</i> OX-19, <i>Morganella morganii</i> IIDKono, <i>Serratia marcescens</i> IFO12648, <i>Enterobacter cloacae</i> IFO13535, <i>Enterobacter aerogenes</i> NCTC10006, <i>Pseudomonas aeruginosa</i> E-2, <i>Acinetobacter calcoaceticus</i> IFO12552.....	Mueller-Hinton agar, 35°C, 24 h
<i>Listeria monocytogenes</i> clinical isolate no. 1, <i>Corynebacterium</i> sp. clinical isolate no. 1.....	Trypticase soy agar with 5% sheep blood, 35°C, 24 h
<i>Streptococcus pyogenes</i> GTC262, <i>Streptococcus agalactiae</i> GTC1234.....	Trypticase soy agar with 5% sheep blood, 35°C, 24 h, 5% CO ₂
<i>Neisseria gonorrhoeae</i> clinical isolate no. 1.....	Thayer-Martin selective agar, 35°C, 24 h, 5% CO ₂
<i>Bacteroides fragilis</i> clinical isolate no. 4.....	Gifu anaerobic medium agar, 35°C, 2 days, anaerobically
<i>Lactobacillus casei</i> ATCC 393, <i>Lactobacillus gasseri</i> JCM1017, <i>Lactobacillus crispatus</i> JCM1030, <i>Lactobacillus acidophilus</i> JCM1132.....	de Man, Rogosa, and Sharpe agar, 35°C, 2 days, anaerobically
<i>Chlamydia trachomatis</i> serovar D/UW-3/Cx.....	Propagation by McCoy cell culture (6)

buffer (see the supplemental material). The extracts were subjected to immunochromatography. The results of these experiments revealed that the minimum cell densities detectable by this method were 10^4 , 10^4 , 10^4 , 10^4 , and 10^5 CFU/ml of *C. albicans*, *C. tropicalis*, *C. guilliermondii*, *C. parapsilosis*, *C. glabrata*, and *C. krusei*, respectively (Fig. 1B [only results obtained with *C. albicans* are shown]). The amount of mannan extracted from 10^4 CFU/ml of *C. albicans* was equivalent to 3.2 ng/ml, as estimated from the above-described calibration curve. Therefore, the amount of mannan in 10^3 CFU/ml should be 0.32 ng/ml, which is below the detection limit of this test, yielding a negative result.

To test the specificity of this immunochromatographic method, bacterial cells at a density of 10^8 CFU/ml were heated and then tested as described above. None of the 28 bacterial species (Table 1) tested yielded a positive result. Similarly, eight fungal species listed in Table 1 were tested with extracts from 10^8 CFU/ml. Among these, only the extract from *Saccharomyces cerevisiae* gave positive results at cell densities over 10^5 CFU/ml but was negative at 10^4 CFU/ml. It was reported that anti-*Candida* antisera cross-reacted with *S. cerevisiae* (1).

To test the reliability of this immunochromatographic test, a total of 200 clinical vaginal swabs were examined and the results were compared with those obtained by the culture method. Of these 200 swabs, 50 were immunochromatography positive; of these 50, 49 were culture positive and 1 was culture negative. The culture method detected *Candida* in swabs from 61 subjects, of which 49 were immunochromatography positive and 12 were negative. Of the 200 subjects, 138 were negative by both immunochromatography and the culture method. Therefore, the total numbers of negative results obtained by immunochromatography and the culture method were 150 and 139, respectively. Immunochromatography yielded false-negative results with 12 samples whose viable *Candida* cell levels were below the limit of detection by the immunochromatography method. One sample showed a very faint line that was barely detectable macroscopically and was regarded as a false-positive result. The reason for this false-positive result is not known. Compared with the results of the culture method, the sensitivity, specificity, positive predictive value, and negative predictive value of the immunochromatography method were calcu-

lated to be 80.3 (49/61), 99.3 (138/139), 98.0 (49/50), and 92.0% (138/150), respectively.

To the best of our knowledge, there has been only one report on an immunochromatographic test for vaginal *Candida* diagnosis (2). Unfortunately, that paper reported only on the testing of clinical subjects without presenting basic study data. The sensitivity, specificity, positive predictive value, and negative predictive value of our method are superior to the values obtained by that method.

Newly developed immunochromatography can detect as few as 10^4 CFU/ml of *Candida* spp. in about 30 min without using expensive equipment or skill. Therefore, the method can be used even in small clinics, as well as large laboratories.

This work was supported by grants-in-aid from the Food Safety Commission Japan.

REFERENCES

1. **Ataoglu, H., J. Zueco, and R. Sentandreu.** 1993. Characterization of epitopes recognized by *Candida* factor 1 and 9 antisera by use of *Saccharomyces cerevisiae* *mnv* mutants. *Infect. Immun.* **61**:3313–3317.
2. **Chatwani, A. J., R. Mehta, S. Hassan, S. Rahimi, S. Jeronis, and V. Dandolu.** 2007. Rapid testing for vaginal yeast detection: a prospective study. *Am. J. Obstet. Gynecol.* **169**:309.e1–309.e4.
3. **Dubois, M., K. A. Gilles, J. K. Hamilton, P. A. Rebers, and F. Smith.** 1956. Colorimetric method for detection of sugars and related substances. *Anal. Chem.* **28**:350–356.
4. **Eckert, L. O.** 2006. Acute vulvovaginitis. *N. Engl. J. Med.* **355**:1244–1252.
5. **Eckert, L. O., S. E. Hawes, C. E. Stevens, L. A. Koutsky, D. A. Eschenbach, and K. K. Holmes.** 1998. Vulvovaginal candidiasis: clinical manifestations, risk factors, management algorithm. *Obstet. Gynecol.* **92**:757–765.
6. **Ikeda-Dantsuji, Y., I. Konomi, and A. Nagayama.** 2005. *In vitro* assessment of the APTIMA Combo 2 assay for the detection of *Chlamydia trachomatis* using highly purified elementary bodies. *J. Med. Microbiol.* **54**:357–360.
7. **Kent, H. L.** 1991. Epidemiology of vaginitis. *Am. J. Obstet. Gynecol.* **165**:1168–1176.
8. **Linhares, I. M., S. S. Witkin, S. D. Miranda, A. M. Fonseca, J. A. Pinotti, and W. J. Ledger.** 2001. Differentiation between women with vulvovaginal symptoms who are positive or negative for *Candida* species by culture. *Infect. Dis. Obstet. Gynecol.* **9**:221–225.
9. **Reiss, E., L. de Repentigny, R. J. Kuykendall, A. W. Carter, R. Galindo, P. Auger, S. L. Bragg, and L. Kaufman.** 1986. Monoclonal antibodies against *Candida tropicalis* mannan: antigen detection by enzyme immunoassay and immunofluorescence. *J. Clin. Microbiol.* **24**:796–802.
10. **Sobel, J. D.** 1997. Vaginitis. *N. Engl. J. Med.* **337**:1896–1903.
11. **Vermitsky, J. P., M. J. Self, S. G. Chadwick, J. P. Trama, M. E. Adelson, E. Mordechai, and S. E. Gyax.** 2008. Survey of vaginal-flora *Candida* species isolates from women of different age groups by use of species-specific PCR detection. *J. Clin. Microbiol.* **46**:1501–1503.